



**Universitätsklinikum
Jena**

Klinik für Hautkrankheiten

Qualitäts-zertifiziert nach DIN EN ISO 9001:2008

Universitätsklinikum Jena · Klinik für Hautkrankheiten · Erfurter Str. 35 · 07743 Jena

smartfiber AG

Anita Varga

Breitscheidstraße 154

D-07407 Rudolstadt

Erfurter Straße 35
07743 Jena

Direktor: Prof. Dr. med. P. Elsner

**Labor für In-vitro-Forschung und
Routinediagnostik**

Laborleiter:

PD Dr. Uta-Christina Hipler

Telefon 03641 93 73 31

Telefax 03641 93 74 37

Ansprechpartner:

Dr. Cornelia Wiegand

Telefon 03641 93 75 84

Telefax 03641 93 74 37

web: [http:// www.derma.uni-jena.de](http://www.derma.uni-jena.de)

25. August 2014

study report:

***In vitro* evaluation of the
anti-oxidative potential of the
textile sample “smartcel™ sensitive”
(# 14 07 088)**

Bachstraße 18 · 07743 Jena · Telefon 03641 93 00

Internet: www.uniklinikum-jena.de

Gerichtsstand Jena

Steuernummer 161 / 144 / 02978 · USt-IdNr. DE 150545777

Bankverbindung:

Sparkasse Jena · BLZ 830 530 30 · Konto 221

Universitätsklinikum Jena · Körperschaft des öffentlichen Rechts

als Teilkörperschaft der Friedrich-Schiller-Universität Jena

Verwaltungsratsvorsitzender: Prof. Dr. Thomas Deufel

Medizinischer Vorstand und Sprecher des Klinikumsvorstandes:

Prof. Dr. Klaus Höffken

Wissenschaftlicher Vorstand: Prof. Dr. Klaus Benndorf

Kaufmännischer Vorstand: Dr. Brunhilde Seidel-Kwem

Index

	page
1. Quality certificate	3
2. General information	
2.1 Test items	4
2.2 Reference items	4
2.3 Sponsor	4
2.4 Test facility	4
2.5 Operating schedule	4
3. GLP and quality assurance statement	5
4. Summary	5
5. Background	6
6. Description of materials and test methods	
6.1 Sample preparation	7
6.2 Determination of the antioxidant capacity against ROS	7
6.1 Determination of the antioxidant capacity against RNS	8
6.1 Statistics	8
7. Deviations from the study protocol	9
8. Archiving	9
9. Results and discussion	10
10. Appendix	
10.1 Abbreviations	11
10.2 Tables and Figures	12
10.3 References	13
10.4 Measurement data	14



1. Quality certificate



Certificate of Registration

QUALITY MANAGEMENT SYSTEM – DIN EN ISO 9001: 2008

This is to certify that

**Universitätsklinikum Jena
Klinik für Hautkrankheiten**
Erfurter Straße 35
07740 Jena
Deutschland

Holds Certificate No: **FS 519135/5409D**

and operates a Quality Management System which complies with the requirements of DIN EN ISO 9001:2008 for the following scope:

Dermatology, Allergology, Occupational Dermatology, Andrology,
Dermatohistology, Dermatooncology, Skin Physiology/Skin Study Center,
Research Laboratory, Diagnostic Laboratory, Dermatological Surgery,
Laser, Photodermatology, Proctology, Phlebology, Wound Healing,
Trichology, Venerology, Pediatric Dermatology

For and on behalf of BSI:



Executive Board, BSI Group Deutschland GmbH

Originally registered: **12.12.2003**

Latest issue: **10.10.2012**

Expiry date: **16.10.2015**

Page: 1 of 1



This certificate was issued electronically and remains the property of BSI and is bound by the conditions of contract. An electronic certificate can be authenticated online. Printed copies can be validated at www.bsi-global.com/ClientDirectory or telephone +49 (0)69 222289-200. Information and Contact: BSI (EMEA) Headquarters: 389 Chiswick High Road, London, W4 4AL, United Kingdom A Member of the BSI Group of Companies.



CONFIDENTIAL – VERTRAULICH

Study report

In vitro evaluation of the anti-oxidative potential of the textile sample "smartcel™ sensitive" (# 14 07 088)



2. General information

2.1 Test items

textile sample "smartcel™ sensitive" (# 14 07 088)

2.2 Reference items

-

2.3 Sponsor

smartfiber AG
Breitscheidstraße 154
D-07407 Rudolstadt
Germany
Person responsible: Anita Varga

2.4 Test facility

Klinik für Hautkrankheiten
Universitätsklinikum Jena
Erfurter Straße 35
D-07740 Jena
Germany
Study director: PD Dr. Uta-Christina Hipler

2.5 Operating schedule

Start of experiments: 2014-08-18
End of experiments: 2014-08-22
Date of final report: 2014-08-25

CONFIDENTIAL – VERTRAULICH

Study report

In vitro evaluation of the anti-oxidative potential of the textile sample "smartcel™ sensitive" (# 14 07 088)



3. GLP and quality assurance statement

I assure that the Test facility complies with the Principles of Good Laboratory Practice. Appropriate and technically valid Standard Operating Procedures are established for the described tests. The Test facility is certified according to DIN EN ISO 9001:2008.

PD Dr. U.-C. Hipler

Laborleiterin

Klinik für Hautkrankheiten

Erfurter Str. 35, D-07740 Jena

Tel.: 03641 / 9-37355

Christina.Hipler@med.uni-jena.de

25.08.2014

Date

Study director: PD Dr. Uta-Christina Hipler

4. Summary

The purpose of this study was to investigate the antioxidative capacity against free radicals such as ROS (reactive oxygen species) and RNS (reactive nitrogen species) *in vitro*.

The textile sample "smartcel™ sensitive" (# 14 07 088) exhibited a significant antioxidative capacity *in vitro*. It could be shown that it is equally effective against ROS and RNS.



5. Background

The skin is the major interface between body and environment. It is the most versatile human organ and plays a key role in protecting the body against environmental influences and participates in the regulation of homeostasis, metabolic processes as well as immunological reactions. Oxidative stress by free radicals accelerates skin aging and has been implicated in dermatological diseases such as atopic dermatitis [Sezer et al. 2007, Briganti & Picardo 2003]. UV light induces the generation of free radicals in the cells; hence, the application of topical antioxidants has been recommended [Masaki 2010, Maela-Azulay & Bagatin 2009]. Textiles are the tissues with the longest contact to the human skin. In the clinically and cosmetically field they are used manifold; one scope of application is for instance the functionalization of fabrics with antioxidants.

The antioxidant capacity (AOC) of soluble substances and other materials, e.g. fabrics, can be monitored and quantified using *in vitro* tests. The several methods are based on different reaction mechanisms and employ various radicals and substrates. Peroxyl radicals (ROO[•]) are the most often used radicals for *in vitro* procedures as they present the key radical for auto oxidation of lipids [Ou et al. 2001]. The ROS (reactive oxygen species) test determines the inhibition of the Pholasin[®] oxidation by superoxide anions and other oxygen radicals. In contrast, the RNS (reactive nitrogen species) test measures the efficacy of an antioxidant to decrease the Pholasin[®] oxidation by peroxynitrite (ONOO⁻). Pholasin[®] is a photo protein isolated from the mollusc *Pholas dactylus*, which emits light in the presence of certain oxidants (chemiluminescence). Previously these tests has been successfully used to determine the antioxidant capacity of wound dressings [Wiegand et al. 2009, Wiegand et al. 2006, Schönfelder et al. 2005] and textiles [Hipler & Wiegand 2011, Wiegand et al. 2010, Fluhr et al. 2010].



6. Description of materials and test methods

6.1 Sample preparation

Textile samples were cut using 8 mm and 5 mm punch biopsies (Stiefel Laboratorium GmbH, Germany) corresponding to 0.5 cm² and 0.25 cm², respectively, and transferred to white 96-well plates (greiner bio-one, Germany).

6.1 Determination of the antioxidant capacity against ROS

The capability of the wound dressings to scavenge free radicals such as ROS (reactive oxygen species) was assessed using the chemiluminescent ABEL[®] Antioxidant Test Kits specific for superoxide anion and other radicals containing Pholasin[®] (Knight Scientific Limited, U.K.).

To each sample the respective assay solutions were added. In brief, 25 µL assay buffer (Lot. GA251A B1 120809), 50 µL Pholasin[®] solution (Lot. AA170A B1 110428) as well as 100 µL solution A (Lot. JA511A A1 130422) were added. Then 25 µL of solution B (Lot. KA611A B3 110912) were injected to each well immediately before measurement. A control without sample was run with each assay. The measurement of luminescence was carried out using the LUMIstar Galaxy plate reader (BMG Labtech GmbH, Germany).

The antioxidant capacity of a sample is expressed as percent reduction of peak luminescence as follows:

$$\% \text{ inhibition} = \frac{[(\text{peak-control}) - (\text{peak-sample})] \times 100}{(\text{peak-control})}$$



6.1 Determination of the antioxidant capacity against RNS

The capability of the wound dressings to scavenge free radicals such as RNS (reactive nitrogen species) was assessed using the chemiluminescent ABEL[®] Antioxidant Test Kits specific for peroxynitrite anion containing Pholasin[®] (Knight Scientific Limited, U.K.).

To each sample the respective assay solutions were added. In brief, 100 µL assay buffer (Lot. TA 508A B2 120416) as well as 50 µL Pholasin[®] solution (Lot. AA 170A A2 130520) were added. Then 50 µL of SIN-1 solution (Lot. UB712B 120330) were injected to each well prior to the measurement. A control without sample was run with each assay. The measurement of luminescence was carried out using the LUMIstar Galaxy plate reader (BMG Labtech GmbH, Germany).

The antioxidant capacity of a sample is expressed as percent reduction of peak luminescence as follows:

$$\% \text{ inhibition} = \frac{[(\text{peak-control}) - (\text{peak-sample})] \times 100}{(\text{peak-control})}$$

6.4 Statistics

Experiments were performed in duplicate and measurements were performed in triplicate. All values are expressed as means ± SD (standard deviation). One-way analysis of variance was carried out to determine statistical significances (Microsoft[®] Excel 2000). Differences are considered statistically significant at a level of $p < 0.05$. Asterisks indicate significant deviations from the control (* $p < 0.05$; ** $p < 0.01$; *** $p < 0.001$).



7. Deviations from the study protocol

There were no deviations from the study protocol.

8. Archiving

The following records will be stored in the archives of the Klinik für Hautkrankheiten, Universitätsklinikum Jena according to the GLP regulations:

A copy of the final report, the study plan and a documentation of all raw data generated during the conduct of the study will be stored for at least 4 years after completion of the study.

Unused test items and reference items are stored for at least 12 month after completion of the study.

Materials and samples that are unstable may be disposed of before that time and without sponsor's prior consent.

Records and reports of the maintenance and calibration of apparatus, validation documentation for computerized systems and the historical file of all Standard Operating Procedures (SOPs) is stored in accordance with the appropriate authorities.



9. Results and discussion

Textile sample “smartcel™ sensitive” (# 14 07 088) exhibited a significant, concentration-dependent capacity to inhibit the formation of free reactive oxygen species (Figure 1). Furthermore, textile sample “smartcel™ sensitive” (# 14 07 088) was able to significantly reduce the formation of reactive nitrogen species (Figure 2).

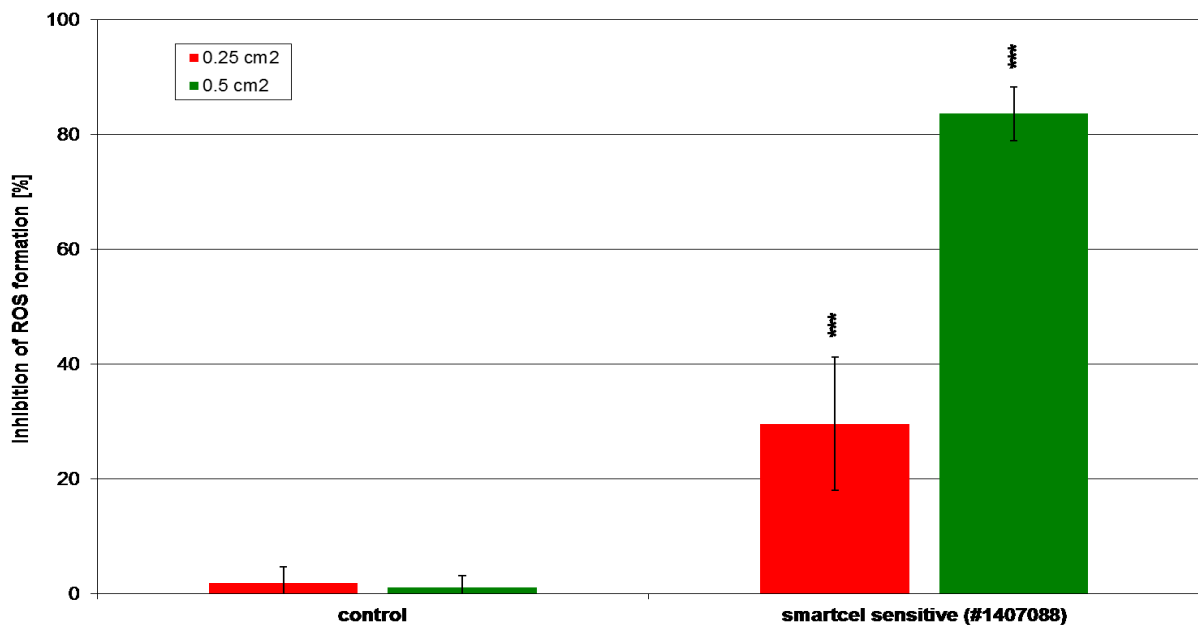


Figure 1: Inhibition of ROS formation (for data see table 1 in the appendix)

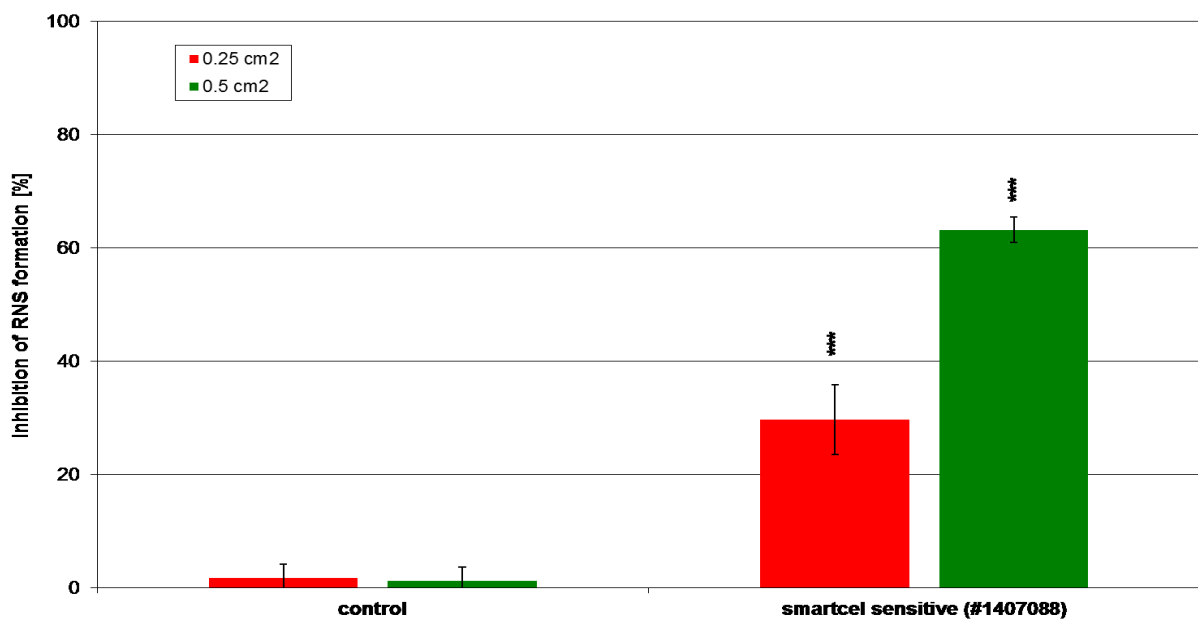


Figure 2: Inhibition of RNS formation (for data see table 2 in the appendix)

CONFIDENTIAL – VERTRAULICH

Study report

In vitro evaluation of the anti-oxidative potential of the textile sample “smartcel™ sensitive” (# 14 07 088)



10. Appendix

10.1 Abbreviations

cm ²	square centimeters
M	molar
mL	milliliters
μL	microliters
mm	millimeters
mM	millimolar
n.s.	not significant
RNS	reactive nitrogen species
ROS	reactive oxygen species
SD	standard deviation



10.2 Tables and Figures

	page
Figure 1: Inhibition of ROS formation (for data see table 1 in the appendix)	10
Figure 2: Inhibition of RNS formation (for data see table 2 in the appendix)	10
Table 1: Reduction of ROS formation.	14
Table 2: Reduction of RNS formation.	14



10.3 References

Briganti S, Picardo M. JEADV 2003; 17:663-9

Fluhr JW, Breternitz M, Kowatzki D, Bauer A, Bossert J, Elsner P, Hipler UC. Exp Dermatol 2010; 19:e9-15

Manela-Azulay M, Bagatin E. Clinics Dermatol 2009; 27:469-474

Masaki H. J Dermatol Sci 2010; 58:85-90

Ou BX, Hampsch-Woodill M, Prior RL. J Agricultur Food Chem 2001; 49:4619-26

Schönfelder U, Abel M, Wiegand C, Klemm D, Elsner P, Hipler U-C. Biomaterials 2005; 26: 6664-6673

Sezer E, Ozugurlu F, Ozyurt H, Sahin S, Etikan I. Clin Exp Dermatol 2007; 32:430-4

Wiegand C, Elsner P, Hipler UC, Klemm D. Cellulose 2006; 13:689-696

Wiegand C, Heinze T, Hipler UC. Wound Rep Reg 2009; 17:511-521

Wiegand C, Fluhr JW, Elsner P, Hipler UC. in Cellulose: Structure and properties, derivatives and industrial uses, eds. A. Lejeune & T. Deprez, Nova Science Publishers, 2010



10.4 Measurement data

Table 1: Reduction of ROS formation.

	sample size	0.25 cm ²	0.5 cm ²
control	[%]	0.00	0.00
		0.00	0.00
		4.98	4.98
		5.78	1.37
		0.00	0.00
		0.00	0.29
	mean	1.8	1.1
	SD	2.8	2.0
textile sample "smartcel™ sensitive" (# 14 07 088)	[%]	42.15	88.62
		30.67	87.95
		17.36	85.73
		28.30	80.37
		15.98	82.43
		42.97	76.46
	mean	29.6	83.6
	SD	11.6	4.7
	p-value (control)	0.0002	0.0001

Table 2: Reduction of RNS formation.

	sample size	0.25 cm ²	0.5 cm ²
control	[%]	0.00	0.00
		0.00	0.00
		6.53	6.53
		0.00	0.00
		0.68	0.00
		2.88	0.90
	mean	1.7	1.2
	SD	2.4	2.4
textile sample "smartcel™ sensitive" (# 14 07 088)	[%]	23.76	62.92
		24.02	62.05
		23.66	65.84
		31.97	61.07
		35.49	66.46
		38.92	60.73
	mean	29.6	63.2
	SD	6.2	2.2
	p-value (control)	0.0001	0.0001